


# Crystalmark Etching Tool - DONATED to UCSB Physics

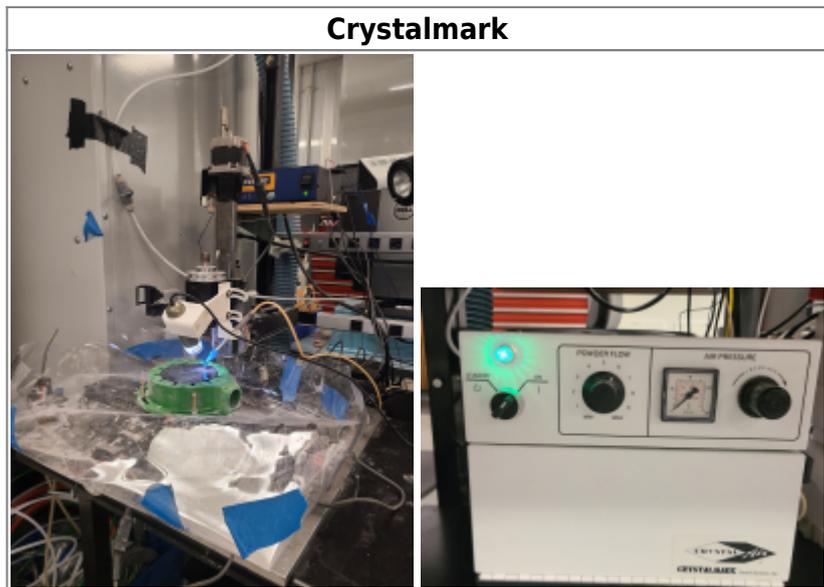

**Tool Type:** CNC Airjet Abrasion Cutter

**Location:** Broida?

**Description:** This is a dental tool meant for cavity prep, which has been repurposed for CNC cutting of glass and silicon. The airjet end of the CrystalMark has been attached to the head of a Sherline CNC mill.

**Manufacturer:** CrystalMark

## About (Legacy - tool disassembled after loss of performance)

The CrystalMark Etcher is located in the Elings 3430 across from the microfluidics station in an acrylic covering. It is attached to the [Sherline CNC Diamond Drill](#).

It uses air pressure and aluminum oxide abrasive to cut or etch patterns in materials such as glass or silicon. The kerf is about 800 microns wide.

Patterns and pre-programmed holes can be uploaded via .dxf files to a program which converts dxf patterns to g code, and holes can also be added manually. (See SOP for detailed instructions)


## Safety Concerns

Inhaling large amounts of aluminum oxide can be hazardous. It is recommended to wear a mask and eye protection when using the CrystalMark.

# Training Documentation

[CrystalMark SOP](#)

## Example Cuts

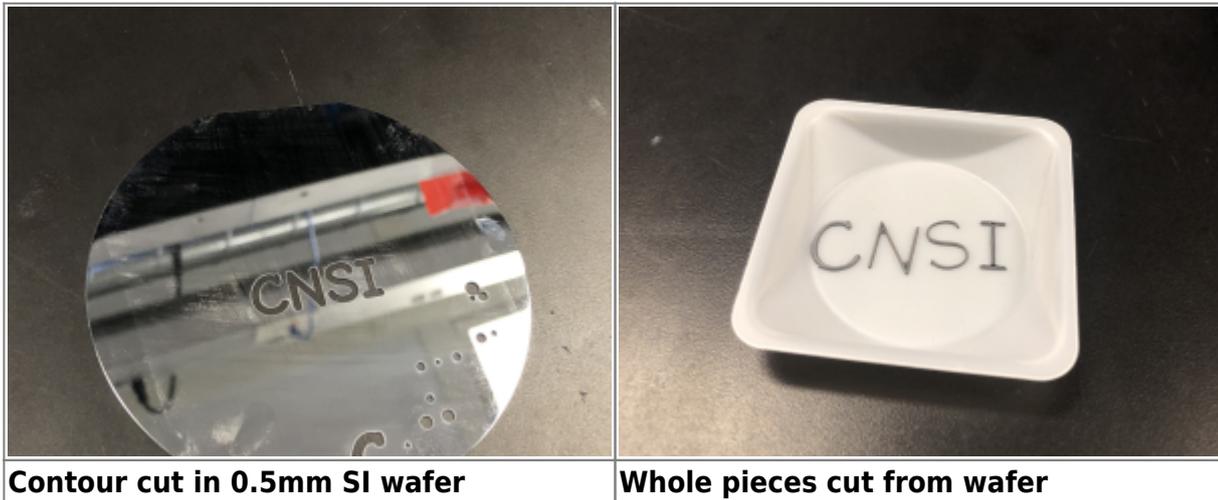

## Detailed Specifications

- Stepper motor mounts and couplers on X-, Y- and Z-axes
- Maximum CNC travel positioning speed: 22 in/min
- Max clearance (table to spindle): 8.00" (203 mm)
- Travel axes x,y,z: 8.65" (220 mm), 5.00" (127 mm), 6.25" (159 mm)
- CNC Stepper motor holding torque: 136 oz-in

## Reference Documentation

Work in progress

From: <https://microfluidics.cnsi.ucsb.edu/wiki/> - Innovation Workshop Wiki

Permanent link: <https://microfluidics.cnsi.ucsb.edu/wiki/doku.php?id=crystalmark>

Last update: **2026/03/17 16:55**



