
Crystalmark Etching Tool - RETIRED

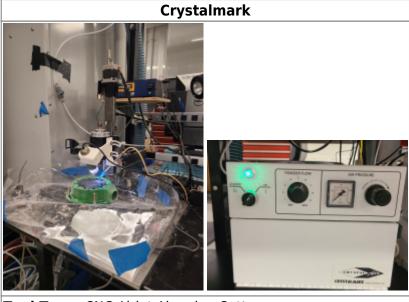

Tool Type: CNC Airjet Abrasion Cutter

Location: Elings 3430

Description: This is a dental tool meant for cavity prep, which has been repurposed for CNC cutting of glass and silicon. The airjet end of the CrystalMark has been attached to the head of a Sherline

CNC mill.

Manufacturer: CrystalMark

About (Legacy - tool disassembled after loss of performance)

The CrystalMark Etcher is located in the Elings 3430 across from the microfluidics station in an arcylic covering. It is attached to the Sherline CNC Diamond Drill.

It uses air pressure and aluminum oxide abrasive to cut or etch patterns in materials such as glass or silicon. The kerf is about 800 microns wide.

Patterns and pre-programmed holes can be uploaded via .dxf files to a program which converts dxf patterns to g code, and holes can also be added manually. (See SOP for detailed instructions)

Safety Concerns


Inhaling large amounts of aluminum oxide can be hazardous. It is recommended to wear a mask and eye protection when using the CrystalMark.

Last update: 2025/06/11 17:03

Training Documentation

CrystalMark SOP

Example Cuts

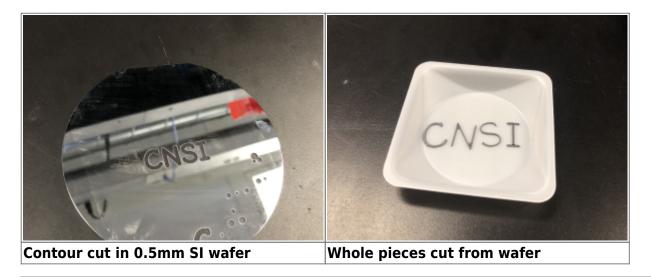

Detailed Specifications

- Stepper motor mounts and couplers on X-, Y- and Z-axes
- Maximum CNC travel positioning speed: 22 in/min
- Max clearance (table to spindle): 8.00" (203 mm)
- Travel axes x,y,z: 8.65" (220 mm), 5.00" (127 mm), 6.25" (159 mm)
- CNC Stepper motor holding torque: 136 oz-in

Reference Documentation

Work in progress

From:

https://microfluidics.cnsi.ucsb.edu/wiki/ - Innovation Workshop Wiki

Permanent link:

https://microfluidics.cnsi.ucsb.edu/wiki/doku.php?id=crystalmark

Last update: 2025/06/11 17:03

