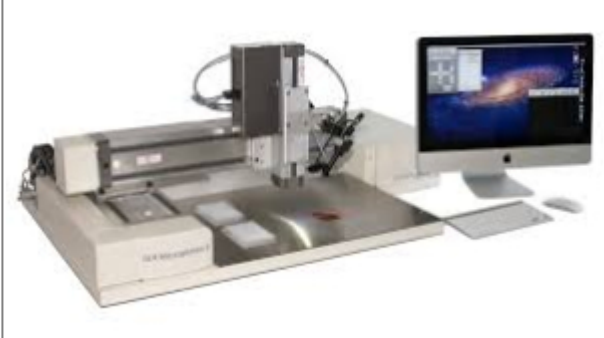


Sonoplot Multimaterial Writer

Sonoplot	
	
Tool Type:	Microplotter
Location:	Microfluidics Lab
Description:	Multimaterial writer
Manufacturer:	Sonoplot

About

The Sonoplot Multimaterial Writer is a microplotting instrument that uses ultrasonic vibrations in its nozzle to deposit small liquid volumes onto planar surfaces in discrete spots or continuous traces. The writer is capable of creating features as small as 5 microns and may be used to create biological microarrays or printed microelectronics.

Safety Concerns

- As with any automated machinery, make sure that your body is clear of the moving parts to avoid injury.
- The glass needle used to fill the liquid well is easily shattered, exercise caution when handling.

Training Documentation

Detailed Specifications

Build Volume: 31x30x9.6 cm

XYZ Resolution: 5 μm

Deposition Volume: ≥ 0.6 pL

Viscosity: ≤ 450 cP

Reference Documentation

MicroPlotter Manual

Micropipette tip replacement guide

From: <https://microfluidics.cnsi.ucsb.edu/wiki/> - Innovation Workshop Wiki

Permanent link: https://microfluidics.cnsi.ucsb.edu/wiki/doku.php?id=sonoplot_multimaterial_writer&rev=1689180004

Last update: **2023/07/12 16:40**

